


Allows under counter access Folding Carboy Spill Stand



Dynalon[®]
Labware

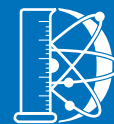

Hinged design to allow bottom 1/2 gallon trough to fold up and away for under counter access. Sliding shelf with thumb holes is removable and accommodates up to 15" diameter carboys.

Custom sizes available upon request.

175 Humboldt Street, Suite 300
Rochester, NY 14610
Telephone 800-334-7585
Fax 585-334-0241
www.dynalon.com

Dynalon Labware is your source for safety products!



Spill Stand allows under counter access due to a fold away spill trough!



Folding Carboy Spill Stand is designed for bench-top use with a folding feature to Allow access to drawers and other materials below. Manufactured of 3/8" High Density Polyethylene, this heavy duty spill tray will hold carboys up to 15" in diameter.

Main spill containment tray has a 2 gallon capacity and bottom spill trough will hold up to 1/2 gallon.


Carboy or other dispensing product sits upon a 16" x 12" sliding tray that moves forward for dispensing and then slides back to allow bottom spill trough to fold up and out of the way.

Spill containment is possible no matter where the tray is positioned.

Custom sizes available upon request.

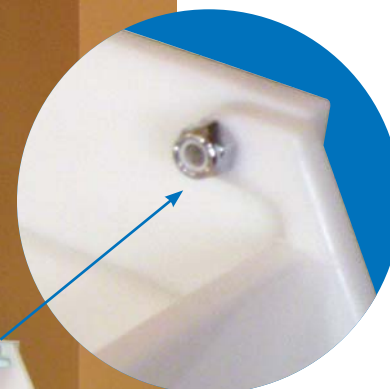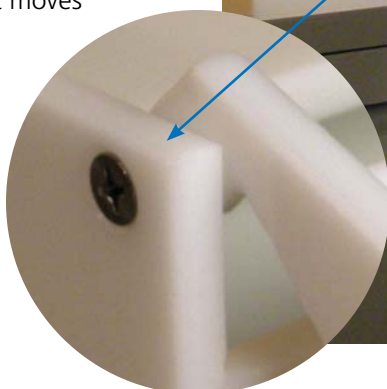

107374	Open	Folded
Depth inches (cm)	21 (53.34)	16 (40.64)
Width inches (mm)	16 (40.64)	16 (40.64)
Height inches (mm)	13.5 (34.29)	10 (25.4)
Weight pounds (Kg)	19 (8.618)	

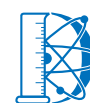

**Dynalon
Labware**

175 Humboldt Street, Suite 300
Rochester, NY 14610
Telephone 800-334-7585
Fax 585-334-0241
www.dynalon.com