


MultiCan® Barrier Containers EVOH Permeation Protection



Dynalon
Labware

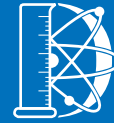

5L



10L



20L

- UN 3H1/Y1 rated
- Ideal for food packaging
- Protection against permeation
- No emission of aromas or smells
- CoEx 6: Co-extrusion up to 6 layers
- Highly effective barrier against gases
- Supplied with regular & Tamper Evident closure

Six-layer multilayer technology MultiCan® Containers





Six-layer multilayer technology makes this container a perfect fit for oxygen-sensitive contents and liquids tending to permeate, such as non-polar hydro- carbon compounds.

Benefits of MultiCan Barrier Containers:

- The filling medium only has contact with food certified HDPE
- The characteristics of your products remains stable
- The permeation of your liquids will be prevented effectively
- Produced according to BRC Standard
- Flexible & inherently stable transport packaging with UN certification

405404	-0005	-0010	-0020
Capacity gallons (liters)	1.32 (5)	2.64 (10)	5.28 (20)
Brimful Capacity gallons (liters)	1.4 (5.3)	3.06 (11.6)	5.94 (22.5)
Weight +5% / -0% g	300	480	920
UN Marking /USA	3H1/Y1.4/150	3H1/Y1.8/150	3H1/Y1.9/150
Pallet Quantity (each)	351	140	54
Pallet Dimensions inches	48 x 42 x 82	46 x 32 x 89	48 x 44 x 52
A1 - Height inches (mm)	9.02 (229)	12.24 (311)	15.71 (399)
A2 - Stack Height inches (mm)	8.5 (216)	11.81 (300)	15.28 (388)
B - Base Width inches (mm)	7.64 (194)	9.13 (232)	11.42 (290)
C - Base Length inches (mm)	6.46 (164)	7.17(192)	9.65 (245)
D1 - Outlet Diameter inches (mm)	2.4 (61)	2.4 (61)	2.4 (61)
E - Height inches (mm)	5.24 (133)	7.32 (186)	8.23 (209)
Thread Neck	DIN 61	DIN 61	DIN 61
Diameter inner neck inches (mm)	1.90 (48.5)	1.90 (48.5)	1.90 (48.5)

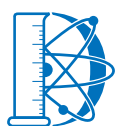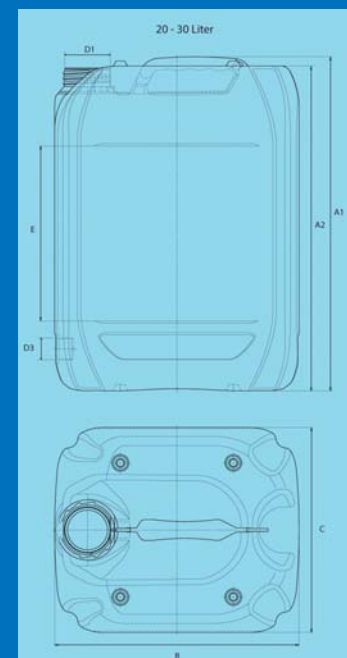

Dynalon® Labware

175 Humboldt St, Ste 300 Rochester, NY 14610
Tel 800.334.7585 www.dynalon.com