

Baritainer® Jerry Cans FDA Compliant & Recyclable







Enhanced chemical resistance Baritainer® Jerry Cans



- Replaces fluorinated, steel and epoxy phenolic
- Enhanced chemical resistance & barrier properties
- Lightweight and stackable
- UN 3H1 Y marking

- Impact resistant
- FDA compliant
- Recyclable

Baritainer® Jerry Cans are the alternative solution for packaging, storing and transporting chemicals and other products that usually require glass, metal, fluorinated HDPE or other single barrier containers. Baritainers can be used to safely transport solvents such as household chemicals, industrial and laboratory chemicals, chemical waste, cleaning solvents, adhesives, wood preservatives, automotive additives, additives in the agricultural, automotive photo chemical sector. They can be used to safely transport foodstuffs such as flavors, edible oils, essential oils and other industry products such as cosmetics and fragrances.

BARITAINER KEY ATTRIBUTES:

KORTRAX® BR/HDPE mixture tight heads can be used to safely transport solvent-based products such as household chemicals, industrial chemicals, cleaning solvents, adhesives, wood preservatives, automotive additives, and agricultural chemicals. They may also be used to safely transport foodstuffs such as flavors, fragrances and essential oils.

BARITAINER SPECIFICATIONS:

Description: UN rated KORTRAX/HDPE mixture extrusion blow-molded tight head with a black tamper evident screw cap

Screw Cap/Neck 4,5 & 10L Baritainers: 50mm LDPE Screw Cap with Tamper

Evident Security Ring and EPDM gasket

Screw Cap/Neck 20L Baritainers: 70mm, 6TPI HDPE Screw Cap, Tamper Evident

Security Ring, with ¾" NPT and EPDM gasket

Pullout Spouts: Available

Colors: White, Natural; custom colors available

BARITAINER COMPOSITION:

Body: KORTRAX® BR and Chemically Resistant HDPE Mixture


Screw Cap: Polyethylene

FDA Statement: All components

are FDA Compliant

HDPE maximum temperature: 248°F (120°C) and minimum temperature: -148°F (-100°C). NOT AUTOCLAVABLE


Baritainer accessories available online.

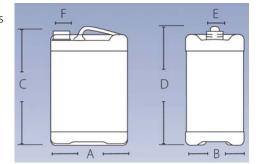

405594	-0004	-0001	-0002	-2070
Capacity gallons (liters)	1.1 (4)	1.32 (5)	2.64 (10)	5.28 (20)
Weight lbs (grams)	0.75-0.84	0.84-0.93	1.68-1.76	2.65
	(340-380)	(380-420)	(760-800)	(1202)
UN Marking /USA	3H1/Y1.8/150	3H1/Y1.8/150	3H1/Y1.2/100	3H1/Y1.8/100
Pallet Quantity (each)	396	360	160	96
A - Length inches (mm)	7.85 (199)	7.85 (199)	9.33 (237)	10.1 (265)
B - Width inches (mm)	6.51 (165)	6.51 (165)	8.43 (214)	10.9 (277)
C - Height inches (mm)	7.13 (181)	8.77 (223)	11.56 (294)	14.3 (363)
E - Cap (mm)	50	50	50	70
F - Neck inches (mm)	1.5 (38)	1.5 (38)	1.5 (38)	2.755 (70)





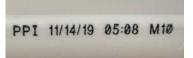



175 Humboldt Street, Suite 300 Rochester, NY 14610 Telephone 800-334-7585 Fax 585-334-0241 www.dynalon.com