



Waste Disposal Containers Bio-bin® and Accessories



Dynalon
Labware



Dynalon Bio-bins® are a cost effective and convenient alternative to rigid plastic waste collection and disposal containers.

Bio-bins are ideal for use in academic, clinical, government & industrial laboratories for a wide range of non-sharps biomedical and laboratory waste that require special heat treatment such as autoclaving or incineration before final disposal.

The ideal solution for biohazard or non-hazardous waste REDUCE Cost, Weight, Transport, CO2 Emissions, Storage Space

- Space efficient, light weight fold-flat design
- Stackable post use. UN 3291 approved
- Printed user instructions & user ID/date information on lid
- Autoclavable. Use with #202065 Autoclave Sterility Tape
- Built-in super absorbent Bio-matt™ to soak up any excess liquid in Bio-bin
- Built-in lid and handles
- Temporary & permanent seals
- Waterproof construction
- Reorder information printed on the Bio-bin
- Supplied with 2 x 2" Biohazard Waste Labels

Catalog Number	Description	UOM
797303-0002	Bio-Bin, Loop Model, 2.5L	50/CS
797303-0006	Bio-Bin, Pipette Model, 6L	40/CS
797303-0030	Bio-Bin, Floor Model, 30L	10/CS
797303-6100	Bio-Bin, U-Stand, for 6L	1 EA
797303-6200	Bio-Bin, Z-Stand, for 6L	1 EA
797303-6300	Bio-Bin, M-Stand, for 30L	1 EA

30 Liter Bio-bin® Spill Containment Trays, PP

- Ideal spill trays for autoclaving waste disposal containers such as the Bio-bins
- Autoclave tray for other plastic or glass laboratory items
- 3/16" double-welded construction
- Help contain accidental chemical spills in the laboratory
- Tested in the USA



Catalog Number	Description	Volume	Dimensions (L x W x H)
797303-6310	Holds one 30L Bio-bin upright	2 Gallon	12 x 12 x 3"
797303-6320	Side Flaps Closed (image above)	3 Gallon	12 x 18 x 3"
797303-6330	Side Flaps Open	4 Gallon	4 x 12 x 26"



Common items such as microbiological loops, pipettes, tips, test tubes, swabs and vials are ideal for disposal using this method.

175 Humboldt Street, Suite 300
Rochester, NY 14610
Telephone 800-334-7585
Fax 585-334-0241
www.dynalon.com

Fold-flat
2.5L and 30L Bio-bin® are delivered folded flat




Temporary Seal
When lid is folded down press here for temporary seal



Permanent Seal
Peel off release paper and fold down lid to seal permanently



Printed user instructions
Simple instructions on how to seal the Bio-bin®

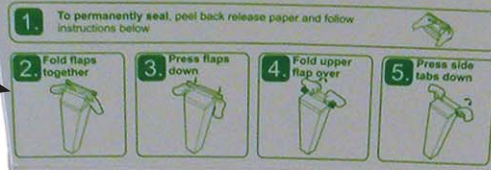

Built in lid
Push inwards here to close lid




Biohazard
Supplied with 2" x 2" Biohazard labels for user application if necessary



Light weight
Weights up to 65% less than plastic alternatives

Water proof construction
96% paper construction with a very thin coating of plastic laminated to make it water proof inside and out. Heat sealed joints prevent any leakage.

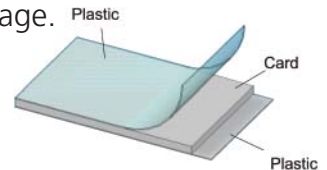


No sharps
Only used for non-sharp infectious waste.



Stackable pre and post use
Once 6L and 30L Bio-bins are sealed down they can be stacked on top of each other.



Absorbent Bio-matt™
Weights up to 65% less than plastic alternatives

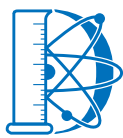

Dynalon Labware