

Customize to suit your needs! Flexible Use Lab Stands



Dynalon[®]
Labware

Funnel
Stand



- Polypropylene chrome-plated steel
- Continuously variable height adjustment with innovative magnet concept
- Adjustable height up to 450 mm

Pipette
Stand




NEW

Sedimentation
Funnel Stand



Variable height adjustment Flexible Use Lab Stands

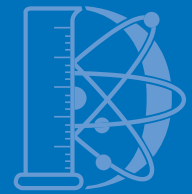

896223-0200 Funnel stand with two funnel holders and two reducing inserts. Space-saving design with two practical storage areas.



896223-0400 Funnel stand with four funnel holders and four reducing inserts. A very sturdy pedestal provides a storage area for four containers and ensures absolute safety against overturning.



871183-2000 Stand with two holders for sedimentation funnels, suitable for funnels 90 mm or more in diameter.



896303-0090 frame for pipettes with two rotating perforated plates and drip pan.

- Polypropylene, chrome-plated steel stands
- Continuously variable height adjustment with innovative magnet concept
- Adjustable height up to 450 mm
- Funnel stands for funnels Ø 40 - 180 mm
- Sedimentation stands for funnels with Ø 90 mm or more
- Pipette stands have 40 slots, Ø 10 - 18 mm
- Additional funnel holders are available as accessories

Sturdy and stable stand, resistant to chemicals, shatter-proof and with continuously variable height adjustment. Thanks to an innovative magnet concept no clamping screws are necessary. The height of the holders can be adjusted easily and flexibly with one hand. The laboratory stand can be flexibly adapted with corresponding accessories, for example with adjustment funnel holders or a holder for Imhoff funnels. The funnel stand can also be converted into a pipette stand.

The **funnel stand** is suitable for all funnels from 40 to 180 mm in diameter. The funnel holders can be used with reducing inserts for small funnels. The funnel stand also makes an excellent filtering stand. 9622-0200 includes frame and two funnel holders with reducing insert, for funnels Ø 40-180 mm. 9622-0400 is supplied with four funnel holders.

The **stand for sedimentation funnels** is used to keep, hold or store Imhoff funnels safely and reliably during sedimentation. Two funnel holders ensure that Imhoff funnels (Ø 40 - 180 mm) will be held safely and securely and can be placed in an absolutely vertical position. 7118-2000 includes frame and two sets of sedimentation funnel holders for Imhoff funnels from Ø 90 mm.

Frame for pipettes with two rotating perforated plates and drip pan. The holders are suitable for up to 40 pipettes of different sizes, for all measuring and full pipettes up to a max. diameter of 18 mm. The holders are suitable for up to 40 pipettes of different sizes, for all measuring and full pipettes up to a max. diameter of 18 mm. A drip pan protects the tips of the pipettes and collects liquid as it drips off. The pipette stand is suitable for drying and for space saving storage of pipettes or thermometers, test tubes and similar items. A very sturdy pedestal ensures absolute safety against overturning. 9630-0090 includes frame and two perforated plates/one drip pan.

Part Number	Description
896223-0200	Funnel stand for two funnels
896223-0400	Funnel stand for four funnels
871183-2000	Imhoff sedimentation funnel stand
896303-0090	Pipette stand
896223-0001	Funnel holder with reducing insert, for funnels Ø 40-180 mm
871183-2001	Sedimentation funnel holder for Imhoff funnels from Ø 90mm
340295	Imhoff Cone, SAN 1000mL
896303-0100	Replacement set for pipette stand, 2 perf. plates/1 drip pan

