




Laboratory Organizer Spill Stand Carboy Spill Containment



Dynalon[®]
Labware

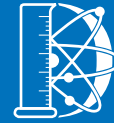

Designed for laboratory carboys and other large containers with spigots to contain liquid spills.

**Custom sizes
available
upon request.**



175 Humboldt Street, Suite 300
Rochester, NY 14610
Telephone 800-334-7585
Fax 585-334-0241
www.dynalon.com



Dynalon Labware is your source for safety products!



Built-in lower bin for convenient storage of labware or small accessories.

- Trays are fabricated with 1/4" thick chemically resistant HDPE. All welded construction featuring a 1/2" skid resistant top shelf
- Provides a safer way to dispense chemicals and other liquids in the lab
- Protects cabinets, laboratory benches and floors from liquid spills
- Each tray contains two convenient drain ports and 1/2" tubulation connection on lower trough
- Built-in side handles for convenient transport
- 8" shelf height to accommodate up to a 5000mL beaker
- All trays are fabricated and tested in our Rochester, New York facility
- **Custom sizes available upon request.**



107364 Carboy Spill Tray Specifications:

Height: 13" / 330.2 mm

Width: 24.5" / 622.3 mm

Depth: 25" / 635 mm

Shelf: 24"W x 12"D / 609 x 305 mm

Tray (lower trough) capacity: 2.5 Gal / 10 L

Weight: 16 lbs / 7.25 Kg

Dynalon®
Labware

quality . value . service

175 Humboldt Street, Suite 300

Rochester, NY 14610

Telephone 800-334-7585

Fax 585-334-0241

www.dynalon.com